



products

USA Manufacturer | Vials • Inserts • Closures • Racks

Choosing a **USA-based manufacturer** for **laboratory vials** offers a range of **critical advantages**

Lab Supplies:


Chromatography, Headspace, and Dram Vials Limited
Volume Inserts, Vial Racks Snap Caps, Plugs, Poly
Crimp™ Seals, Closures, Environmental Sampling
Bottles and Jars, VOA Vials, Deep Well Storage and
Collection Plates



Big Mouth Crimp Top Vials, 12x32mm, 11mm Crimp Finish

Finneran Products Certified For Science™

Vials

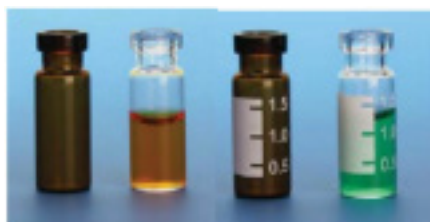

320111-1232A, 320111-1232,
32111E-1232A, 320111E-1232

40% larger opening prevents broken needles due to increased target area. Neck finish allows use of patented Poly Crimp™ Seals or standard aluminum seals. Designed to work in Agilent and other autosamplers. Choose from clear or amber Type I borosilicate glass. ID design incorporates the Step Vial feature that precisely centers a limited volume insert in the vial neck. Available with graduated marking spots.

Shelf Pack - 100 pieces; Case Pack - 1000 pieces

Cat. No.	Description
32011L-1232	2.0mL Clear Vial
32011L-1232A	2.0mL Amber Vial
32011L-1232Z	2.0mL Clear Vial, Silanized
32011L-12AZ	2.0mL Amber Vial, Silanized
32011LE-1232	2.0mL Clear Vial, White Graduation Spot
32011LE-1232A	2.0 Amber Vial, White Graduation Spot

Limited Volume Big Mouth Crimp Top Vials

Tapered Bottom Vials - Manufactured from Type I borosilicate clear glass. The tapered bottom allows for the concentration of samples into the tip by solvent evaporation.

Case Pack - 100 pieces



31111L-1232

Cat. No.	Description
31111L-1232	1.1mL Clear Tapered Bottom Vial

Interlocked Vials

Next Generation Interlocked Vial is made of clear borosilicate glass and features an enhanced insert that provides better sample extraction. Sample capacity of 300µL, 12mm in diameter, 32mm in height, and 11mm large opening crimp neck finish. Vial compatible with all our 11mm seals or closures.

Shelf Pack - 100 pieces; Case Pack - 1000 pieces



30211LNG-12, 30211LMNG-12,
3011LNG-12A

Cat. No.	Description
30211LNG-12	300uL Clear Vial
30211LMNG-12	300uL Clear Vial, with Marking Spot
30211LNG-12A	300uL Amber Vial

Fox Scientific, Inc.

8221 East FM 917 • Alvarado, TX 76009 • Phone:(800)369-5524 • Fax:(817)783-3571

www.foxscientific.com • sales@foxscientific.com

Big Mouth Crimp Top Vials, 12x32mm, 11mm Crimp Finish

Closures

Using Patented Poly Crimp™ seals, in place of aluminum seals, eliminates the need for crimpers and decappers

Poly Crimp™ Seals - Patented

The Finneran patented polypropylene Poly Crimp™ seal fits virtually any 12x32mm sample vial with a crimp style finish. The Poly Crimp™ seal snaps on and off with crimping being an option for a tighter seal. Poly Crimp™ seals are available in a variety of colors for easy sample identification.

Shelf Pack - 100 pieces; Case Pack - 1000 pieces

Cat. No.	Description
5500-11*	11mm Clear Poly Crimp™ Seal, w/Molded Septum
5500SB-11*	11mm Clear Poly Crimp™ Seal, w/Starburst
5510-11*	11mm Clear Poly Crimp™ Seal, w/10mil PTFE
5517-11*	11mm Clear Poly Crimp™ Seal, w/10mil PTFE w/Slit
5540-11*	11mm Clear Poly Crimp™ Seal, Clear PTFE/Butyl Rubber
5550-11*	11mm Clear Poly Crimp™ Seal, PTFE/Silicone
5560-11*	11mm Clear Poly Crimp™ Seal, PTFE/Silicone/PTFE
5570-11*	11mm Clear Poly Crimp™ Seal, PTFE/Silicone w/Slit
5581-11*	11mm Clear Poly Crimp™ Seal, PTFE/Silicone w/Starburst

* Available in additional colors; (B) Blue, (G) Green, (R) Red, & (Y) Yellow



Aluminum Seals, Lined

Made of the highest quality aluminum; for use with standard 11mm chromatography vials. The open top allows for filling and retrieving samples with a syringe. Seals are available in a variety of colors for easy sample identification.

Shelf Pack - 100 pieces; Case Pack - 1000 pieces

Cat. No.	Description
5110-11*	11mm Silver Seal, 10mil PTFE
5120-11	11mm Silver Seal, 10mil Polypropylene
5140-11*	11mm Silver Seal, Clear PTFE/Red Rubber
5140-11X*	11mm Silver Seal, Clear PTFE/Natural Red Rubber
5141-11	11mm Silver Seal, Clear PTFE/Butyl Rubber
5150-11*	11mm Silver Seal, PTFE/Silicone
5160-11*	11mm Silver Seal, PTFE/Silicone/PTFE
5170-11	11mm Silver Seal, PTFE/Silicone w/Slit
5190-11	11mm Silver Seal, Black Viton® Lined

* Available in additional colors; (B) Blue, (G) Green, (R) Red, & (Y) Gold



11mm Magnetic Seals

Magnetic seals for use with instruments (CTC Analytics PAL) transporting the vials using magnets.

Shelf Pack - 100 pieces; Case Pack - 1000 pieces

Cat. No.	Description
5140M-11	11mm Silver Magnetic Seal, PTFE/Red Rubber
5150M-11	11mm Silver Magnetic Seal, PTFE/Red Silicone
5160M-11	11mm Silver Magnetic Seal, PTFE/Silicone/PTFE



Fox Scientific, Inc.

8221 East FM 917 • Alvarado, TX 76009 • Phone: (800) 369-5524 • Fax: (817) 783-3571

www.foxscientific.com • sales@foxscientific.com

J.G. Finneran Associates, Inc. — Visit www.JGFinneran.com • 1-800-552-3696 • 1-856-696-3605

Engineered and manufactured in compliance with ISO 9001:2015 certified guidelines